



The Levanto is a safe and comfortable standing and transfer device designed to help patients with some trunk balance and weakened leg muscles to stand up with a natural motion pattern. It features clip sling attachments and its compact design makes it the ideal option for smaller spaces.

#### Optimal Features

The Levanto is an electric sit-to-stand device that encourages a natural standing motion and is designed for standing, sitting, dressing, toileting and balance training for users with low core or leg strength. The shin pad angle adjusts with the patient while they are being raised for optimal comfort and security and the legs are easily spread electronically to allow direct access while seated (for example, in a wheelchair). The non-slip surface footplate offers stable support and is removable for easier cleaning and hygiene control.

#### Ideal for Home Care Use

The Levanto features a compact footprint that allows it to easily pass through narrow

doorways. It features a lifting arm with additional handles for comfort and easy usage. The shorter legs allow for easier rotation, especially in a home environment that may have less space for manoeuvring. With a LINAK battery and actuator, you are ensured years of trouble-free use.

#### Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

#### HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B  
117 43 Stockholm  
Phone:- +46 8 665 35 00  
Fax:- +46 8 665 35 10  
info.se@humancaregroup.com

#### HUMAN CARE CANADA

10-155 Colonnade Road  
Ottawa, ON K2E 7K1  
Phone: +613.723.6734  
Fax: +613.723.1058  
info.ca@humancaregroup.com

#### HUMAN CARE NEDERLAND

Elspeterweg 124  
8076 PA Vierhouten  
Phone: +31 577 412 171  
Fax: +31 577 412 170  
info.nl@humancaregroup.com

#### HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K  
Austin, TX 78754  
Phone: +512.476.7199  
Fax: +512.476.7190  
info.us@humancaregroup.com

## Accessories



### Slings

Assortment of sizes and designs to accommodate most every lifting need.



### Battery

Complete system with control box enclosed in a discreet module. Included with lift.



### Hand Control

Slow or fast setting for arm up/down, electrical leg spreading. Included with lift.



### Wheels

Rear wheels feature kick-brake. Diameter 125mm.

## Lifting Situations



### Use to/from a bed

Use to raise/lower a patient from/to a bed.



### Use to/from a chair

Use to raise/lower a patient from/to a chair.



### Use to/from the toilet

Use to raise/lower a patient from/to a toilet.

## Features




### Control Box

Emergency stop, emergency lowering and soft-start functions.




### Lifting Arm

Lift arm extends 80cm / 31.5" to accommodate taller patients.

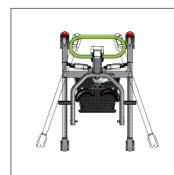

### Spreadable Base

Electrical adjustment.

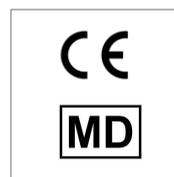


Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

## Technical Specifications

Levanto			
Part No.	20000170	IP Rating	IPX4
Width	68.7cm / 27"	Turning Radius	128cm / 50.3"
Length	110cm / 43.3"	Front Wheels	10cm/4"
		Rear Wheels with Brake	12.5cm/5"
Height	Arm Lowered: 88.9cm / 35" Arm Raised <sup>1</sup> : 164.8cm / 64.9"	Battery	24V DC/2.9Ah, 2x 12V/2.9Ah
Height from Floor	11.9cm / 4.7"	Weight	49kg / 108lb
Adjustable Leg (Min - Max)	Outside: 67.9cm - 113.1cm / 26.7" - 44.5"	Weight Capacity	160kg / 352.7lb
Sound Level	48dB (loaded) while lifting 52dB (loaded) while spreading legs		

<sup>1</sup>Measurement based on arm mast at highest setting

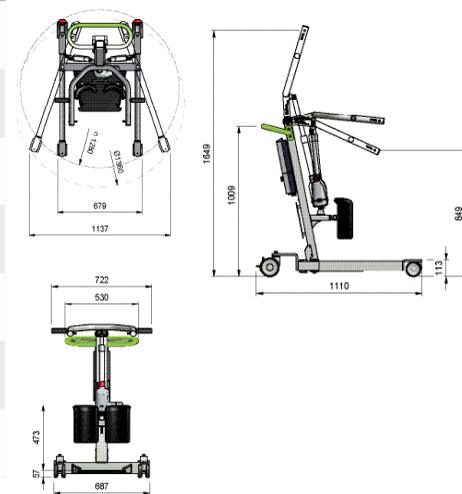

*We improve freedom of movement and quality of care for people.*

[www.humancaregroup.com](http://www.humancaregroup.com)

**HUMAN CARE SWEDEN (HQ)**  
 Årstaängsvägen 21B  
 117 43 Stockholm  
 +46 8 665 35 00  
 info@humancaregroup.com

**HUMAN CARE CANADA**  
 10-155 Colonnade Road  
 ON K2E 7K1 Ottawa  
 +1 613 723 6734  
 info.ca@humancaregroup.com

**HUMAN CARE THE NETHERLANDS**  
 Elspeterweg 124  
 8076 PA Vierhouten  
 +31 577 41 21 71  
 info.nl@humancaregroup.com

**HUMAN CARE UNITED STATES**  
 8006 Cameron rd, Suite K  
 TX 78754 Austin  
 +1 5 124 76 7199  
 info.us@humancaregroup.com